

# WiTouch<sup>®</sup> Pro

BLUETOOTH<sup>®</sup> TENS THERAPY

## Quick Start Guide

Smartphone not required for operation

**INSTRUCTIONS MUST BE FOLLOWED. PLEASE DO NOT DISCARD.**



For **SETUP & Bluetooth<sup>®</sup>**  
**PAIRING INSTRUCTIONS**


Scan QR code with camera on smartphone  
or visit [hollywog.com/setup](http://hollywog.com/setup)

**⚠ IMPORTANT!** Refer to the user manual for all Indications, Contraindications, Warnings, Precautions, Adverse Reactions, and directions. A detailed user manual with instructions for use is available in the WiTouch Pro App's information tab and for download at [hollywog.com/contact](http://hollywog.com/contact). A printed copy of the user manual can be requested by contacting Hollywog at [info@hollywog.com](mailto:info@hollywog.com).

### Box Contents



WiTouch<sup>®</sup> Pro Bluetooth<sup>®</sup> device (1)

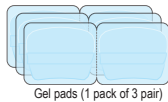

Gel pads (1 pack of 3 pair)



#0 Phillips head  
screwdriver (1)



AAA alkaline  
batteries (2)

## Download, install, & connect the WiTouch Pro App

App is **ONLY** compatible with **iOS 9.0 & above OR ANDROID 4.4 & above.**

- 1 Download the WiTouch Pro App.



Location **must** be **ALLOWED**.\*

- 2 Using the screwdriver provided, remove the battery cover and insert the **AAA** batteries.



- 3 Turn device **ON** by pressing the **ON/OFF** button. The green LED will begin flashing, indicating the device is in standby mode.



- 4 Open your WiTouch Pro app.



- 5 Click **"OK"** for iOS or **"ALLOW"** for Android to enable Bluetooth.



*must be turned on in your  
smartphone settings to pair.*

- 6 Press & HOLD the **+** & **-** intensity buttons at the same time until the LED stops flashing and remains solid green. Let go of the buttons.  
This should take approximately 3-4 seconds **ONLY**.

- 7 Click **CONTINUE** at the bottom of your screen.

- 8 Click **"PAIR"** for iOS device. Click **"OK"** to pair for Android devices.

**NOTE:** If **DEVICE NOT FOUND** repeat steps 6-8.

## Smartphone Pairing

- 1 WiTouch Pro can only be paired with one (1) smartphone at a time, or pairing will be lost.

- 2 **IF YOUR DEVICE WILL STILL NOT PAIR...** ensure Bluetooth is ON, hold both PLUS (+) and MINUS (-) buttons for approximately 10 seconds (hold past the solid green). After approximately 10 seconds your LED will blink Yellow. Let go of the buttons. Bluetooth is now ON and ready to be paired (Refer to steps 6-8 above to follow pairing directions). If you see a green flash instead of yellow, **REPEAT** the approximate 10-second hold. The Yellow blink indicates Bluetooth is ON and ready to be paired.


- 3 Refer to Page 5 in the User Manual at [hollywog.com/contact](http://hollywog.com/contact) or call (423) 305-7777.

**FLIP ME OVER!**



## Apply Gel Pads


1 Remove one (1) pair of gel pads, and separate the two gel pads.

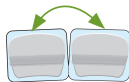

2 Remove & discard the blue two-piece liner.




3 Adhere the exposed side of the gel pad to the black electrode surface on the back of your unit. To ensure good adhesion, rub firmly along the surface to remove large air bubbles.

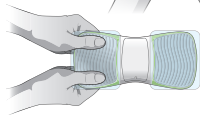


Repeat steps 2 & 3 for the second gel pad.

## Place on Your Back

1 Remove the protective liners from the gel pads.

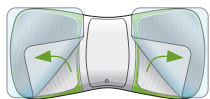

2 Align the center of the device over the spine and press the device to your back.

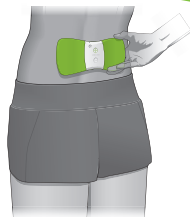

## Start a Treatment

1 Within the App, choose one of the four treatment programs and press START. Your 30-minute Treatment will begin.

2 Adjust the intensity using **+** and **-** controls within the app, or **+** and **-** on the device itself.



## Using the WiTouch Pro WITHOUT a Smartphone

1 Follow the APPLY GEL PADS instructions above.

2 Turn device ON by pressing the ON/OFF button. The green LED will begin flashing, indicating the device is in standby mode.

3 Follow the PLACE ON YOUR BACK instructions above.

4 To start the WiTouch Pro Exclusive Treatment Program press the (+) button. To increase the intensity, continue to press the (+) button one press at a time. Pressing the (-) button (one press at a time) will decrease the intensity.

Further instructions are available on Page 7 of the User Manual @ [hollywog.com/contact](http://hollywog.com/contact).

**Hollywog**  
100 Cherokee Blvd.  
Chattanooga, TN 37405 USA  
[www.hollywog.com](http://www.hollywog.com)

© and © 2021 Hollywog, LLC. All rights reserved. The Bluetooth® word mark and logos are registered trademarks owned by the Bluetooth SIG, Inc. and any use of such marks by Hollywog, LLC is under license. Other trademarks and trade names are those of their respective owners.  
U.S. Patents 9,764,133, 8,972,016, D712,052, D716,958. International patents pending. Made in USA